



Laparoscopy Hand Instruments

It's your choice.



LAPAROSCOPY HAND INSTRUMENTS, 360° ROTATABLE



THE STANDARD The modular system for whatever you need

- Long-life and reliable materials
 Large variety of handle designs
 Ratchet mechanism / deactivation


More than 400 different laparoscopic inserts (Ø 3.0 / 3.5 / 5.0 / 10mm)



The RZ spherical head locking system for fast assembling / disassembling



Modular 3-piece detachable instrument system consists of:

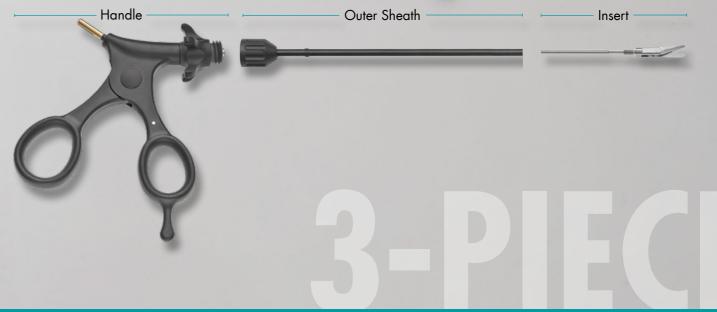

Please find more information in our Laparoscopy catalogue





LAPAROSCOPY HAND INSTRUMENTS, 360° ROTATABLE



THE ENDO TIP SYSTEM

The environmental solution for the always sharp scissors

- Reusable handle and shaft, disposable tip
 Always having a sharp tip
 Ecological and economical reduces waste and saves budget



Complete Instrument consists of:

- PEEK Polymer Handle (300-450-300)
 Actuator rod (300-450-006)
- Outer tube including screw cap blue, \emptyset 5 mm (300-450-005) 8
- Insert

PEEK POLYMER HANDLES



PEEK Polymer Handle, blue rotating wheel 300-450-300 without ratchet

The EndoTip system allows to reuse the handle and the sheath of the Instrument, but gives the opportunity to easily change the tip. That provides an always sharp scissors without creating the waste of a complete single use instrument.

4 METZENBAUM

non-sterile, limited reusable medium, 15mm curved 300-450-177 short, 10mm, curved 300-450-171 long, 18 mm, curved 300-450-173 long, 18 mm, straight 300-450-174 medium, 15mm, straight 300-450-176

METZENBAUM

sterile, single use medium, 15mm curved short, 10mm, curved

300-450-177S 300-450-1715





LAPAROSCOPY HAND INSTRUMENTS, 360° ROTATABLE



SINGLE USE SHEATH SYSTEM Convenient and sharp

- Combine your preferred handle with a sharp singe use scissors
 No need for resharpening of your most preferred Metzenbaum
 High quality handle combined with single use scissors



DISPOSABLE PART

PEEK Polymer Handle, with HF connection 300-750-101 without ratchet

> Please find more information in our Laparoscopy catalogue



The fully disposable sheath system gives you the opportunity to combine your prefered MIC instrument handle with a disposable sheath. That safes space on your tray and reduces time in reprocessing of the equipment.

METZENBAUM

Scissors Inserts, sterile, box with 10 pcs., Ø 5mm, 330mm working length medium, 15mm curved

300-600-177S

METZENBAUM non-steril

straight, curved, s curved, straight, straight, curved,

| le, limited reusable, 1 pc. | |
|-----------------------------|-------------|
| short, 10mm | 300-600-170 |
| hort, 10mm | 300-600-17 |
| 18 mm, long | 300-600-173 |
| 18 mm, long | 300-600-174 |
| 15mm, medium | 300-600-176 |
| 15mm, medium | 300-600-177 |
| | |





RZ Medizintechnik GmbH

Unter Hasslen 20 78532 Tuttlingen, Germany

Phone: +49 (0)7462/9470-0 Fax: +49 (0)7462/9470-50

sales@rz-medizintechnik.com www.rz-medizintechnik.com

