

VERSATIP

The versatile
solution
for Endo-Urology



VERSATIP



The next generation instruments in Endo-Urology

- ▼ Detachable – one handle for all inserts
- ▼ Color code – the easy identification of the tip
- ▼ Two handles – self holding and standard
- ▼ Easy connection – the reliable ball tip technique
- ▼ Light weight handle – the perfect fit for your fingers
- ▼ Make cleaning easy – flushing port for the complete instrument
- ▼ 360° rotation - just with your finger tip
- ▼ All parts can be ordered seperately



Watch the
product video





360°



Multiple tips design



Color coded inserts



Unique, detachable system



Flushing port for easy cleaning



Two handles - multiple options

FORCEPS & PUNCHES



- ▶ Automatic self-holding
- ▶ Black color coded



351-100-550
Self-holding handle for
forceps and punches
(recommended handle)



Grasping Forceps with toothed cup,
double action, 600mm working length

			Insert only	Complete instrument
	semi-rigid	3 Charr., Ø 1.0mm	253-950-503	253-950-003
		4 Charr., Ø 1.3mm	253-950-504	253-950-004
		5 Charr., Ø 1.6mm	253-950-505	253-950-005
	flexible	3 Charr., Ø 1.0mm	253-951-503	253-951-003
		4 Charr., Ø 1.3mm	253-951-504	253-951-004
		5 Charr., Ø 1.6mm	253-951-505	253-951-005



Stone Grasping Forceps for large stones,
double action jaws, 600mm working length

Insert only

Complete instrument



semi-rigid

3 Charr., Ø 1.0mm

253-952-503

253-952-003

4 Charr., Ø 1.3mm

253-952-504

253-952-004

5 Charr., Ø 1.6mm

253-952-505

253-952-005



flexible

3 Charr., Ø 1.0mm

253-953-503

253-953-003

4 Charr., Ø 1.3mm

253-953-504

253-953-004

5 Charr., Ø 1.6mm

253-953-505

253-953-005



HAMADA long Stone Grasping Forceps,
double action jaws, 600mm working length

Insert only

Complete instrument



semi-rigid

3 Charr., Ø 1.0mm

253-954-503

253-954-003

4 Charr., Ø 1.3mm

253-954-504

253-954-004

5 Charr., Ø 1.6mm

253-954-505

253-954-005



flexible

3 Charr., Ø 1.0mm

253-955-503

253-955-003

4 Charr., Ø 1.3mm

253-955-504

253-955-004

5 Charr., Ø 1.6mm

253-955-505

253-955-005



Biopsy Forceps, 600mm working length

Insert only

Complete instrument



semi-rigid

3 Charr., Ø 1.0mm

253-956-503

253-956-003

4 Charr., Ø 1.3mm

253-956-504

253-956-004

5 Charr., Ø 1.6mm

253-956-505

253-956-005



flexible

3 Charr., Ø 1.0mm

253-957-503

253-957-003

4 Charr., Ø 1.3mm

253-957-504

253-957-004


5 Charr., Ø 1.6mm

253-957-505

253-957-005

VERSATIP

SCISSORS

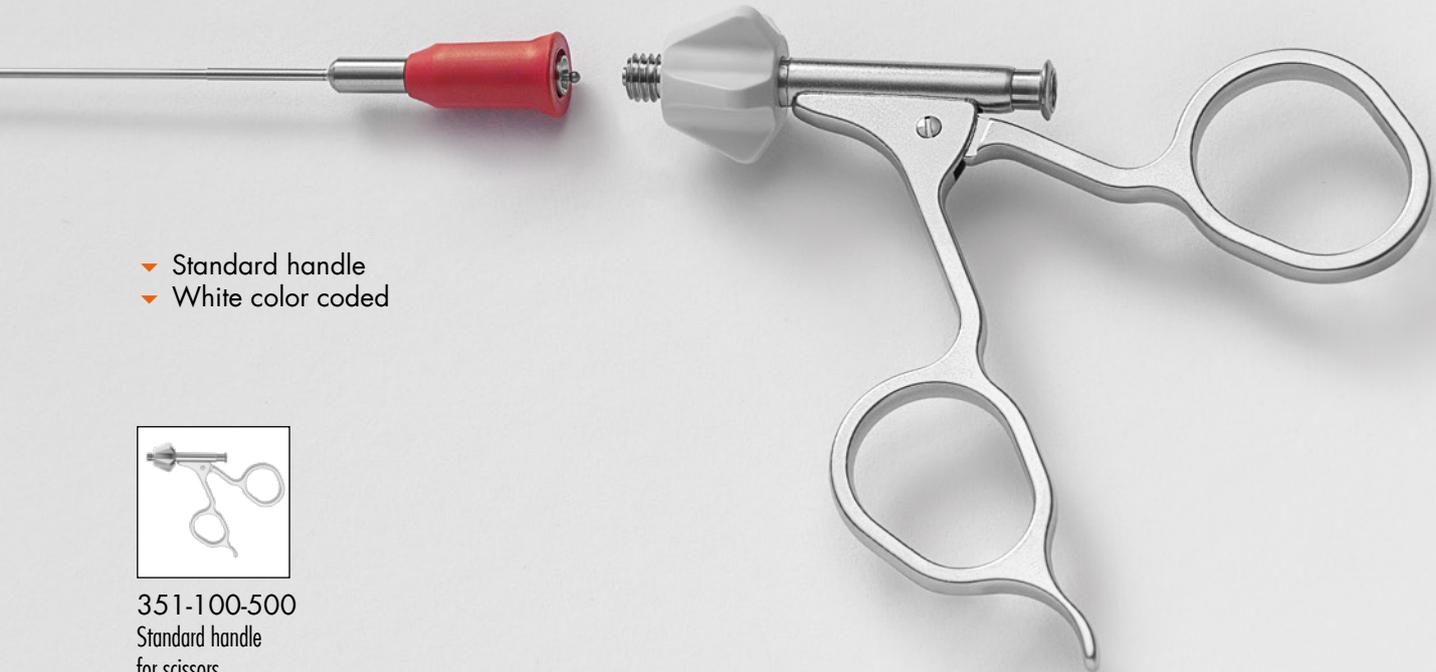

- Standard handle
- White color coded




351-100-500
Standard handle
for scissors
(recommended handle)




Scissors, single action, 600mm working length		Insert only	Complete instrument
 semi-rigid	3 Charr., Ø 1.0mm	253-958-503	253-958-003
	4 Charr., Ø 1.3mm	253-958-504	253-958-004
	5 Charr., Ø 1.6mm	253-958-505	253-958-005
 flexible	3 Charr., Ø 1.0mm	253-959-503	253-959-003
	4 Charr., Ø 1.3mm	253-959-504	253-959-004
	5 Charr., Ø 1.6mm	253-959-505	253-959-005

ASSEMBLY

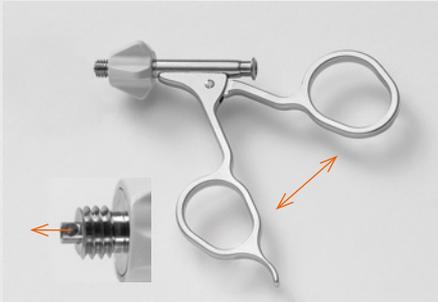

STEP 1

Open the handle.
The ball socket pushes out.
Keep the handle open.

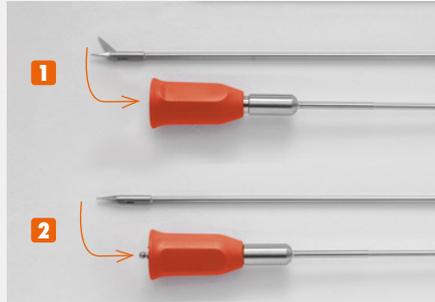

STEP 2

1 Opened jaw
2 Close the jaw - the ball head pushes out. Keep the jaw closed.




STEP 3

Push back the locking nut.



STEP 4

With the handle opened and the jaws kept closed, Insert the ball head into the socket.

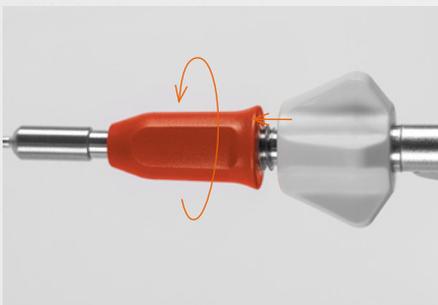

STEP 5

Close the handle. Shaft and insert will be retracted.
Keep the handle closed.
Screw the locking nut onto the handle.



DONE

Check the connection for being tight.
The knob rotates 360° and the handle as well as the jaw opens and closes.



RZ Medizintechnik GmbH

Unter Hasslen 20/22
78532 Tuttlingen, Germany

+49 7462 9470-0

sales@rz-medizintechnik.com
www.rz-medizintechnik.com

